

STAGE 1

Local anaesthesia


Use Anesthetic PRO FORMULA topical anaesthetic. Apply a dense layer to the skin, cover with a film to enhance the anaesthetic effect, exposure time 15–20 min. Remove residual anaesthetic agent with Decasan.

STAGE 2

Setting up the device




1 Remove the Perfoskin Device's portable body (a).

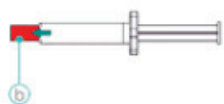


2 Remove the rubber cap (b) from the Perfoskin Prederm syringe.




3 Insert the syringe into Perfoskin Device (c).

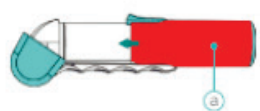

4 Slide on the portable body (a).

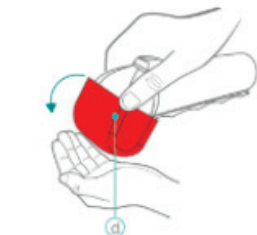

5 Press the edge of the cap (d) with your thumb to remove it. Be careful not to injure yourself.

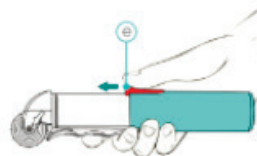

6 Squeeze out a small amount of product by gently pressing the shutter (e).

Perform the procedure by simultaneously pressing the shutter while sliding the roller over the skin surface.

STAGE 3

PERFORMANCE OF PERFOSKIN COMPLEX PROCEDURE


During the first minute of the procedure, apply slight pressure to the roller so that the patient becomes accustomed to the sensation. Then, continue the procedure, increasing roller pressure. (pinpoint bleeding is normal and necessary).

During the procedure, one should cover small areas of skin, respecting the Langer tension lines.

(Fig. 1) Roll the device over the skin in multiple directions based on the problem.

Anti-age, reduced skin turgor

— horizontal and vertical movements.

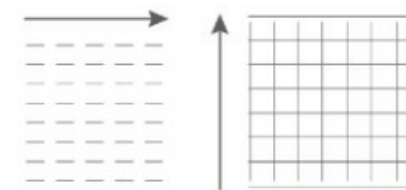

Scars, post-acne, hyperpigmentation

— horizontal, vertical, and diagonally in both directions.

